




Product catalog
CAD / CAM
Materials

CLOSE TO **YOU**

Welcome to **your** Universe

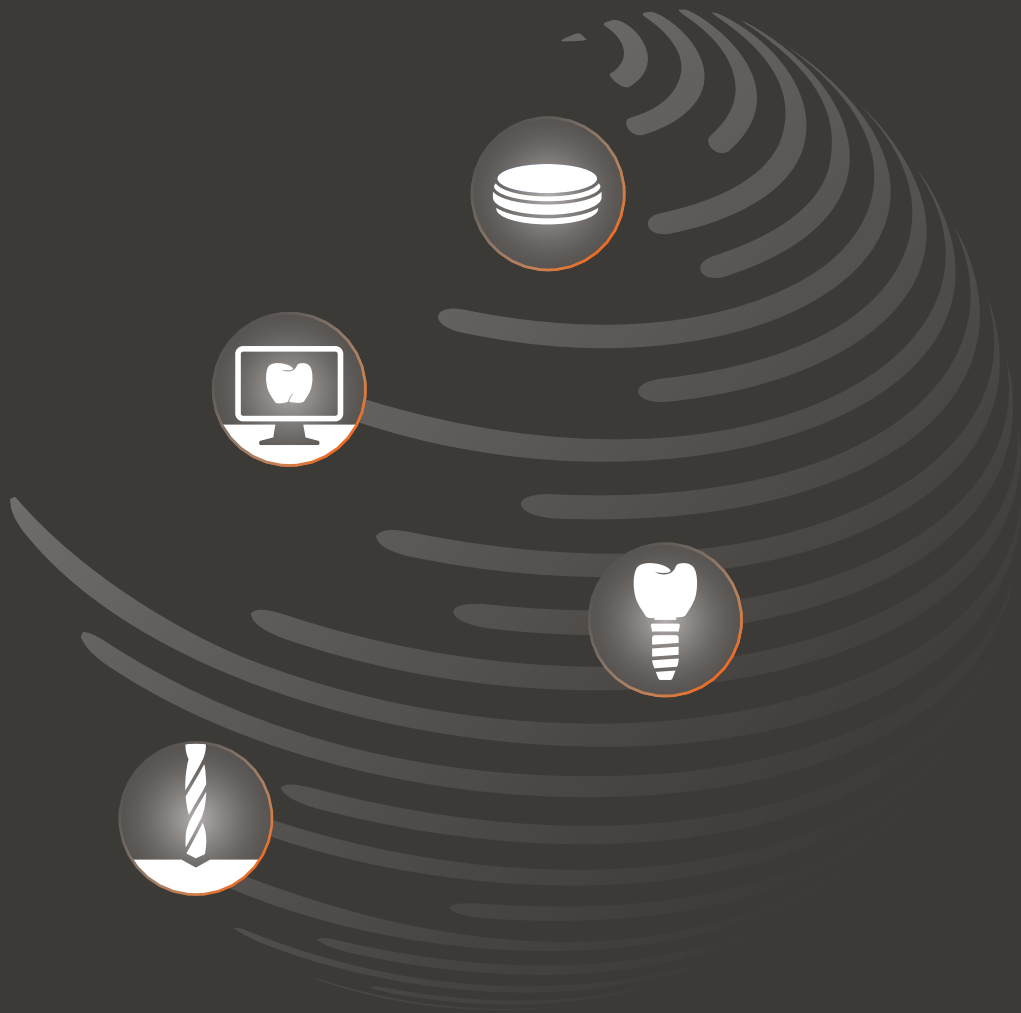

Materials

Technologies

Implant Components

Milling Services

01 High Performance Ceramics		4
	Raw Materials	5
	All zirconia you need	6
	Aesthetics overview	8
	Indication overview	10
Super High Translucent (5Y-TZP)	DD cubeX ² ®	12
Super High Translucent (6Y-PSZ) High Translucent (3Y-TZP)	Nacera® Pearl Natural	16
High Translucent Plus (4Y-TZP)	DD cube ONE®	20
High Translucent (3Y-TZP-LA)	DD Bio ZX ²	26
High Strength (3Y-TZP-A)	DD Bio Z	30
02 Coloring Liquids		35
	Coloring Liquid recommendation	36
	DD Basic Shade	38
	DD Pro Shade	39
	DD Art Elements	41
	DD Shade Concept®-Sets	42
	DD Incisal X	43
03 Aesthetic Finish		44
	DD contrast® color and form concept	45
	DD contrast® set	48
	panther – The Premium Finish	49
04 High Performance Polymers		50
	DD Bio Splint P HI	51
	DD Bio Splint C	52
	DD Bio Splint FLEX	53
	DD provi P HI	54
	DD base P HI	55
	DD poly X ML	56
	DD peek MED	58
05 More Materials		59
	DD Bio CW	60
	DD cam WAX	61
	DD calibration disc	61
06 Accessories		62
	DD Shade Guide	63
	DD Solid Link	64
	DD Solid Link Try-in	65
	DD Uni Bond	66
	DD Uni Bond D	67
	DD phoeniX	68

High Performance Ceramics

01

The raw materials

For us, it is a matter of course that only materials with the best properties come into consideration for dental products. The biocompatibility is the ultimate priority, in the selection of raw materials. Moreover only a high-purity, nanoscale ZrO₂ is suitable for medical use. To ensure reproducible batch stability, the raw material must indicate an extremely homogenous grain size distribution.

Chemical composition [wt.%]

	DD Bio Z (color) 3Y-TZP-A	DD Bio ZX² (color) 3Y-TZP-LA	DD cube ONE[®] (ML) 4Y-TZP	Nacera[®] Pearl Natural 6Y-PSZ + 3Y-TZP	DD cubeX²[®] (ML) 5Y-TZP
ZrO ₂ + HfO ₂ + Y ₂ O ₃	≥ 99,0	≥ 99,0	≥ 99,0	≥ 99,0	≥ 99,0
Y ₂ O ₃	< 6	< 6	< 8	< 10	< 10
Al ₂ O ₃	< 0,5	≤ 0,15	< 0,15	0,10	≤ 0,01
Other oxides	< 1	< 1	< 1	< 1	< 1

	Bio Z color 3Y-TZP-A	DD Bio ZX² color 3Y-TZP-LA	DD cube ONE[®] ML 4Y-TZP	Nacera[®] Pearl Natural 6Y-PSZ + 3Y-TZP	DD cubeX²[®] ML 5Y-TZP
Color method Industrial precolored blanks	Multi Additive Technology [®]	Multi Additive Technology [®]	Multi Additive Technology [®]	–	Multi Additive Technology [®]
Additive	Fe ₂ O ₃ Er ₂ O ₃ MnO ₂	Fe ₂ O ₃ Er ₂ O ₃ MnO ₂	Fe ₂ O ₃ Er ₂ O ₃ Co ₃ O ₄	Fe ₂ O ₃ Er ₂ O ₃ MnO ₂	Fe ₂ O ₃ Er ₂ O ₃ MnO ₂
Colors	in 8 VITA [®] shades + 5 Universal shades	in 16 VITA [®] shades + 5 Universal shades	in 16 VITA [®] shades + 1 Bleach shade	in 16 VITA [®] -Farben + 3 Bleach shades	in 16 VITA [®] shades + 1 Bleach shade





Multi Additive Technology[®] – more than just color







We achieve brilliant color effects by adding selected ions of rare soils or transition metals into the zirconium matrix. Our balanced ion formulation and the homogeneous distribution in the zirconia structure is what we call Multi Additive Technology[®]. It is the key to excellent optical characteristics.

By means of our Multi Additive Technology colored blanks show a natural chroma based on the definition in the VITA classical[®] color scale.










All zirconia you need

	Super High Translucent (5Y-TZP)	Super High Translucent (6Y-PSZ) High Translucent (3Y-TZP)
<p>MULTI LAYER</p> 	<p>DDcubeX²ML</p>  <p>16 VITA®- + 1 Bleach shade</p>	<p>Nacera® Pearl Natural</p>  <p>16 VITA®- + 3 Bleach shades</p>
<p>MONO CHROME</p> 		
<p>WHITE</p> 	<p>DDcubeX²</p> 	
<p>Dentine Coloring Liquids</p>	<p>DDProShade[©] DDBasicShade</p>	

<p>High Translucent Plus (4Y-TZP)</p>	<p>High Translucent (3Y-TZP-LA)</p>	<p>High Strength (3Y-TZP-A)</p>
<p>DDcubeONE^{ML}</p>  <p>16 VITA[®]- + 1 Bleach shade</p>		
	<p>DDBioZX²color</p>  <p>16 VITA[®]- + 5 Universal shades</p>	<p>DDBioZcolor</p>  <p>8 VITA[®]- + 5 Universal shades</p>
<p>DDcubeONE[®]</p> 	<p>DDBioZX²</p> 	<p>DDBioZ</p> 
<p>DDProShade^Z DDBasicShade</p>	<p>DDProShade^Z DDBasicShade</p>	<p>DDProShade^Z DDBasicShade</p>

Aesthetics overview

Generation	Multilayer	
	16 VITA® - + 1 Bleach shade type II, class 4, ≤ 3 units	16 VITA® - + 3 Bleach shade type II, class 5, ≥ 4 units
Strength		
	700 (±100) MPa*	1.000 MPa* – incisal 1.200 MPa* – cervical
Translucency		
	Super High Translucent (5Y-TZP)	Super High Translucent (6Y-PSZ) High Translucent (3Y-TZP)
Masking	slight stump discoloration 	medium stump discoloration
	aesthetic-line yellow-warm	aesthetic-line blue-cold
Aesthetics		
	Processing	automated

16 VITA® - + 1 Bleach shade

type II, class 5, ≥ 4 units



Monochrome

16 VITA® - + 5 Universal shades

type II, class 5, ≥ 4 units



8 VITA® - + 5 Universal shades

type II, class 5, ≥ 4 units



1.200 (±150) MPa*

1.150 (±150) MPa*

1.150 (±200) MPa*



High Translucent Plus (4Y-TZP)



High Translucent (3Y-TZP-LA)



High Strength (3Y-TZP-A)



strong stump discoloration



aesthetic-line yellow-warm



aesthetic-line yellow-warm



aesthetic-line yellow-warm



manual

*measured according to DIN EN ISO 6872

Indication overview

		veneers	inlays/ onlays	reduced crown	monolithic crown anterior tooth
Dental Direkt					
		<i>DDcubeX²</i> <i>DDcubeX²ML</i>			
Nacera Pearl Natural					
<i>DDcubeONE</i> <i>DDcubeONE_{ML}</i>					
<i>DDBioZX²</i> <i>DDBioZX²color</i>					
<i>DDBioZ</i> <i>DDBioZcolor</i>					

*An individual selection of the appropriate material is necessary to ensure the best possible restoration.

Influence on the final shade of a fixed restoration

When choosing a zirconia material for a specific clinical case several aspects should be considered, because the individual patients conditions might offer very different preconditions.

An important factor is the shade of the prepared tooth. If there is no serious deviation from the desired final shade of the fixed restoration, a super high translucent zirconia (DD cubeX²) might be the material of choice.

The quantity of substance removal by preparation can be significant as well, as it affects the wall thickness of the restoration. There is a direct context between the underlying tooth structure, the wall thickness and the resulting shade effect of the restoration.

monolithic CROWN posterior tooth	Hybrid ABUTMENT CROWN	monolithic BRIDGE (up to 3 units)	monolithic BRIDGE	reduced BRIDGE (up to 3 units)	reduced BRIDGE	Hybrid ABUTMENT	BAR

Indication guide_DD medical zirconia_Rev.07_2022/08

In case of a significantly discolored stump and preparation with low substance removal a more opaque zirconia (DD Bio ZX², DD Bio Z) could be used to cover up the underlying structure.


Especially when working with a high translucent plus (DD cube ONE[®]) or super high translucent (DD cubeX²[®]) zirconia material also the cement has to be well-chosen, taking its level of opacity and shade into consideration.

Opaque or colored cementation materials might change the final shade of a restoration.



You find our cementation guide on www.dentaldirekt.de/en/downloads

DDcubeX²[®]

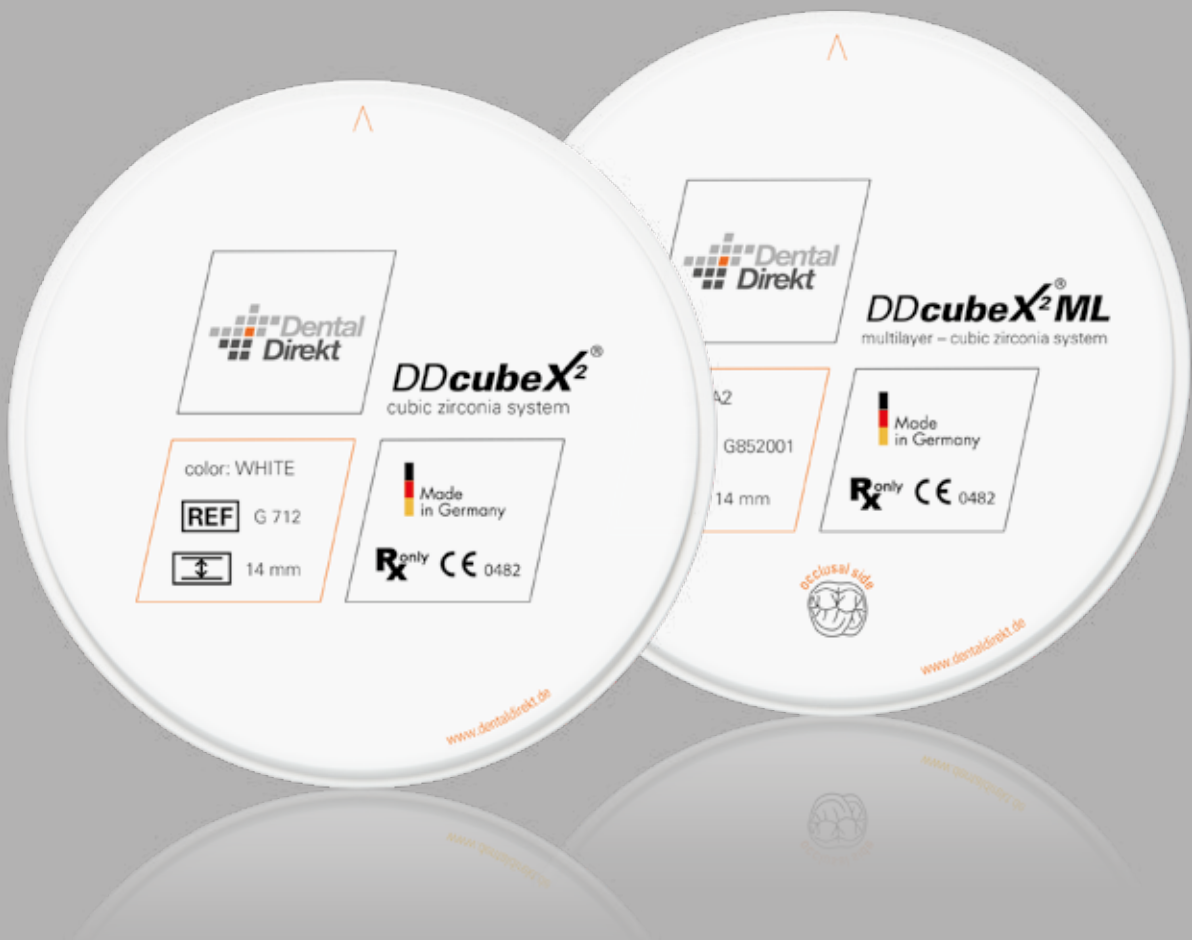

DDcubeX²[®]

cubic zirconia system

Super High Translucent (5Y-TZP)

WHITE

Benefits

- More stable than lithium disilicate with comparable aesthetics
- Our favorite for the anterior region

Perfectly suitable for

- High aesthetic monolithic crowns and bridges (up to 3 units)
- Especially for monolithic anterior restorations (bridges up to 3 units)
- Cut back (one layer)
- High aesthetic veneering



Technical data

Type	5Y-TZP
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 4 – max. 3 unit bridges
Color	White
Aesthetics	Super High Translucent
Sintering temperature	1.450 °C
Flexural strength	700 (±100) MPa*
CAM system	Open systems for Ø 98,5 mm blanks
DD Coloring liquids	DD Pro Shade C DD Basic Shade DD Art Elements
Aesthetic Finish	DD contrast[®]
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872

Measure	Order-no.	Price (euro/net)
98,5 x 10 mm	G 710	164,00
98,5 x 12 mm	G 711	184,00
98,5 x 14 mm	G 712	195,00
98,5 x 16 mm	G 713	215,00
98,5 x 18 mm	G 714	236,00
98,5 x 20 mm	G 715	246,00
98,5 x 25 mm	G 716	287,00



Consult instructions for use.

Not all products are available in every country.

DDcubeX²ML[®]

multilayer – cubic zirconia system

Super High Translucent (5Y-TZP)



Benefits

- More stable than lithium disilicate with comparable aesthetics
- Our favorite for the anterior region in 16 VITA[®]- + 1 Bleach shade
- Incisal layer and dentine layer + 2 transition layers = natural color gradient

Perfectly suitable for

- High aesthetic monolithic crowns and bridges (up to 3 units)
- Especially for monolithic anterior restorations (bridges up to 3 units)
- Cut back (one layer)
- High aesthetic veneering



Technical data

Type	5Y-TZP
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 4 – max. 3 unit bridges
Color	Precolored acc. to VITA [®]
Aesthetics	Super High Translucent
Sintering temperature	1.450 °C
Flexural strength	700 (±100) MPa*
CAM system	Open systems for Ø 98,5 mm blanks
Multi Additive Technology [®]	Reliable reproduction of color via industrial coloring
DD Coloring liquids	DD Art Elements
Aesthetic Finish	DD contrast[®]
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872



Measure	Bleach	A1	A2	A3	A3,5	A4	Price (euro/net)
98,5 x 14 mm	G852017	G852007	G852001	G852002	G852003	G852009	195,00
98,5 x 18 mm	G854017	G854007	G854001	G854002	G854003	G854009	236,00
98,5 x 22 mm	G856017	G856007	G856001	G856002	G856003	G856009	287,00



Consult instructions for use.

Not all products are available in every country.



Measure	B1	B2	B3	B4	Price (euro/net)
98,5 x 14 mm	G852008	G852004	G852010	G852011	195,00
98,5 x 18 mm	G854008	G854004	G854010	G854011	236,00
98,5 x 22 mm	G856008	G856004	G856010	G856011	287,00




Measure	C1	C2	C3	C4	Price (euro/net)
98,5 x 14 mm	G852012	G852005	G852013	G852014	195,00
98,5 x 18 mm	G854012	G854005	G854013	G854014	236,00
98,5 x 22 mm	G856012	G856005	G856013	G856014	287,00



Measure	D2	D3	D4	Price (euro/net)
98,5 x 14 mm	G852015	G852006	G852016	195,00
98,5 x 18 mm	G854015	G854006	G854016	236,00
98,5 x 22 mm	G856015	G856006	G856016	287,00

Nacera® Pearl Natural

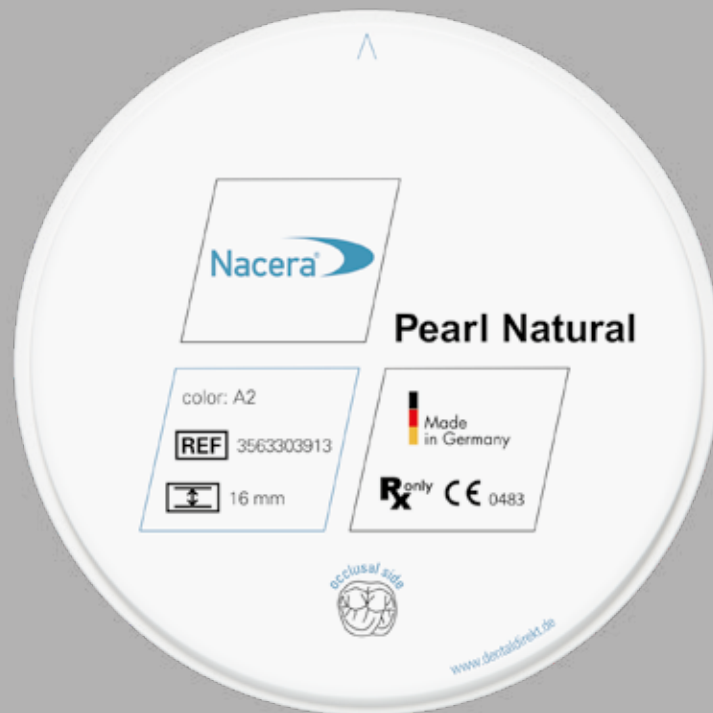

Nacera® Pearl Natural

Super High Translucent (6Y-PSZ)
High Translucent (3Y-TZP)



Benefits

- Naturally soft color and translucency transitions for lifelike looking monolithic restorations (of any span range)
- 16 VITA® shades, two VITA® bleach shades and one special bleach shade
- Highest precision due to distortion-free material
- High process reliability during processing – uncritical positioning during nesting
- Fast and cost-effective for natural-looking monolithic restorations

Perfectly suitable for

- All indications: from monolithic single crowns to complex superstructures for the anterior and posterior region



Technical data

Type	6Y-PSZ + 3Y-TZP
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges with max. 2 connected pontics
Color	Precolored acc. to VITA®
Aesthetics	incisal: Super High Translucent cervical: High Translucent
Sintering temperature	1.500 °C
Flexural strength	incisal: 1000 MPa* cervical: 1200 MPa*
CAM system	Open systems for Ø 98,5 mm blanks
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872



Consult instructions for use.

Not all products are available in every country.

Nacera® Pearl Natural

Super High Translucent (6Y-PSZ)

High Translucent (3Y-TZP)



Measure	Stellar White	OM2	OM3	Price (euro/net)
98,3 x 16 mm	3559303913	3560303913	3561303913	215,00
98,3 x 20 mm	3559303096	3560303096	3561303096	287,00
98,3 x 25 mm	3559303097	3560303097	3561303097	339,00



Measure	A1	A2	A3	A3,5	A4	Price (euro/net)
98,3 x 16 mm	3562303913	3563303913	3564303913	3565303913	3566303913	215,00
98,3 x 20 mm	3562303096	3563303096	3564303096	3565303096	3566303096	287,00
98,3 x 25 mm	3562303097	3563303097	3564303097	3565303097	3566303097	339,00



Consult instructions for use.

Not all products are available in every country.



Measure	B1	B2	B3	B4	Price (euro/net)
98,3 x 16 mm	3567303913	3568303913	3569303913	3570303913	215,00
98,3 x 20 mm	3567303096	3568303096	3569303096	3570303096	287,00
98,3 x 25 mm	3567303097	3568303097	3569303097	3570303097	339,00




Measure	C1	C2	C3	C4	Price (euro/net)
98,3 x 16 mm	3571303913	3572303913	3573303913	3574303913	215,00
98,3 x 20 mm	3571303096	3572303096	3573303096	3574303096	287,00
98,3 x 25 mm	3571303097	3572303097	3573303097	3574303097	339,00



Measure	D2	D3	D4	Price (euro/net)
98,3 x 16 mm	3575303913	3576303913	3577303913	215,00
98,3 x 20 mm	3575303096	3576303096	3577303096	287,00
98,3 x 25 mm	3575303097	3576303097	3577303097	339,00

DDcubeONE®

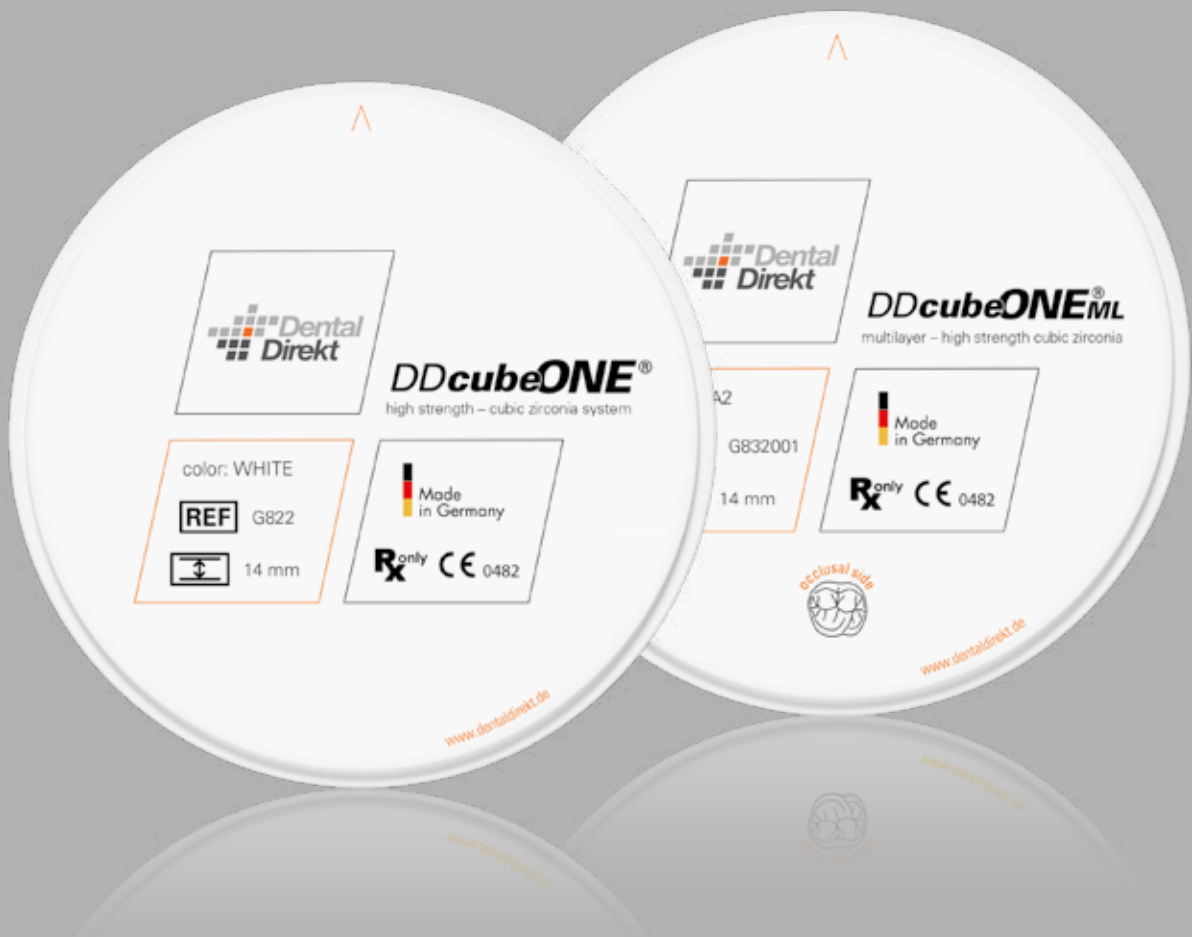

DDcubeONE®

high strength – cubic zirconia system

High Translucent Plus (4Y-TZP)

WHITE

Benefits

- Multi-indicative universal zirconia
- Ideal balance of stability and aesthetics
- Extremely high fracture toughness

Perfectly suitable for

- High aesthetic monolithic crowns and bridges (of any span range)
- Especially for monolithic anterior restorations
- Cut back (one layer)
- High aesthetic veneering



Technical data

Type	4Y-TZP
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges with max. 2 connected pontics
Color	White
Aesthetics	High Translucent Plus
Sintering temperature	1.450 °C
Flexural strength	1.200 (±150) MPa*
CAM system	Open systems for Ø 98,5 mm blanks
DD Coloring liquids	DD Pro Shade Z DD Basic Shade DD Art Elements
Aesthetic Finish	DD contrast®
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872

Measure	Order-no.	Price (euro/net)
98,5 x 10 mm	G820	164,00
98,5 x 12 mm	G821	169,00
98,5 x 14 mm	G822	195,00
98,5 x 16 mm	G823	199,00
98,5 x 18 mm	G824	236,00
98,5 x 20 mm	G825	246,00
98,5 x 25 mm	G826	287,00



Consult instructions for use.

Not all products are available in every country.

DDcubeONE[®]ML

multilayer – high strength cubic zirconia

High Translucent Plus (4Y-TZP)



Benefits

- Multi-indicative universal zirconia in 16 VITA®- + 1 Bleach shade
- Incisal layer and dentine layer + 2 transition layers = natural color gradient
- Ideal balance of stability and aesthetics

Perfectly suitable for

- High aesthetic monolithic crowns and bridges (of any span range)
- Especially for monolithic anterior restorations
- Cut back (one layer)
- High aesthetic veneering



Technical data

Type	4Y-TZP
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges with max. 2 connected pontics
Color	Precolored acc. to VITA®
Aesthetics	High Translucent Plus
Sintering temperature	1.450 °C
Flexural strength	1.200 (±150) MPa*
CAM system	Open systems for Ø 98,5 mm blanks
Multi Additive Technology®	Reliable reproduction of color via industrial coloring
DD Coloring liquids	DD Art Elements DD Incisal X
Aesthetic Finish	DD contrast®
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872



Measure	Bleach	A1	A2	A3	A3,5	A4	Price (euro/net)
98,5 x 14 mm	G832017	G832007	G832001	G832002	G832003	G832009	195,00
98,5 x 18 mm	G834017	G834007	G834001	G834002	G834003	G834009	236,00
98,5 x 22 mm	G836017	G836007	G836001	G836002	G836003	G836009	287,00



Consult instructions for use.

Not all products are available in every country.



Measure	B1	B2	B3	B4	Price (euro/net)
98,5 x 14 mm	G832008	G832004	G832010	G832011	195,00
98,5 x 18 mm	G834008	G834004	G834010	G834011	236,00
98,5 x 22 mm	G836008	G836004	G836010	G836011	287,00




Measure	C1	C2	C3	C4	Price (euro/net)
98,5 x 14 mm	G832012	G832005	G832013	G832014	195,00
98,5 x 18 mm	G834012	G834005	G834013	G834014	236,00
98,5 x 22 mm	G836012	G836005	G836013	G836014	287,00



Measure	D2	D3	D4	Price (euro/net)
98,5 x 14 mm	G832015	G832006	G832016	195,00
98,5 x 18 mm	G834015	G834006	G834016	236,00
98,5 x 22 mm	G836015	G836006	G836016	287,00

DDBioZX²

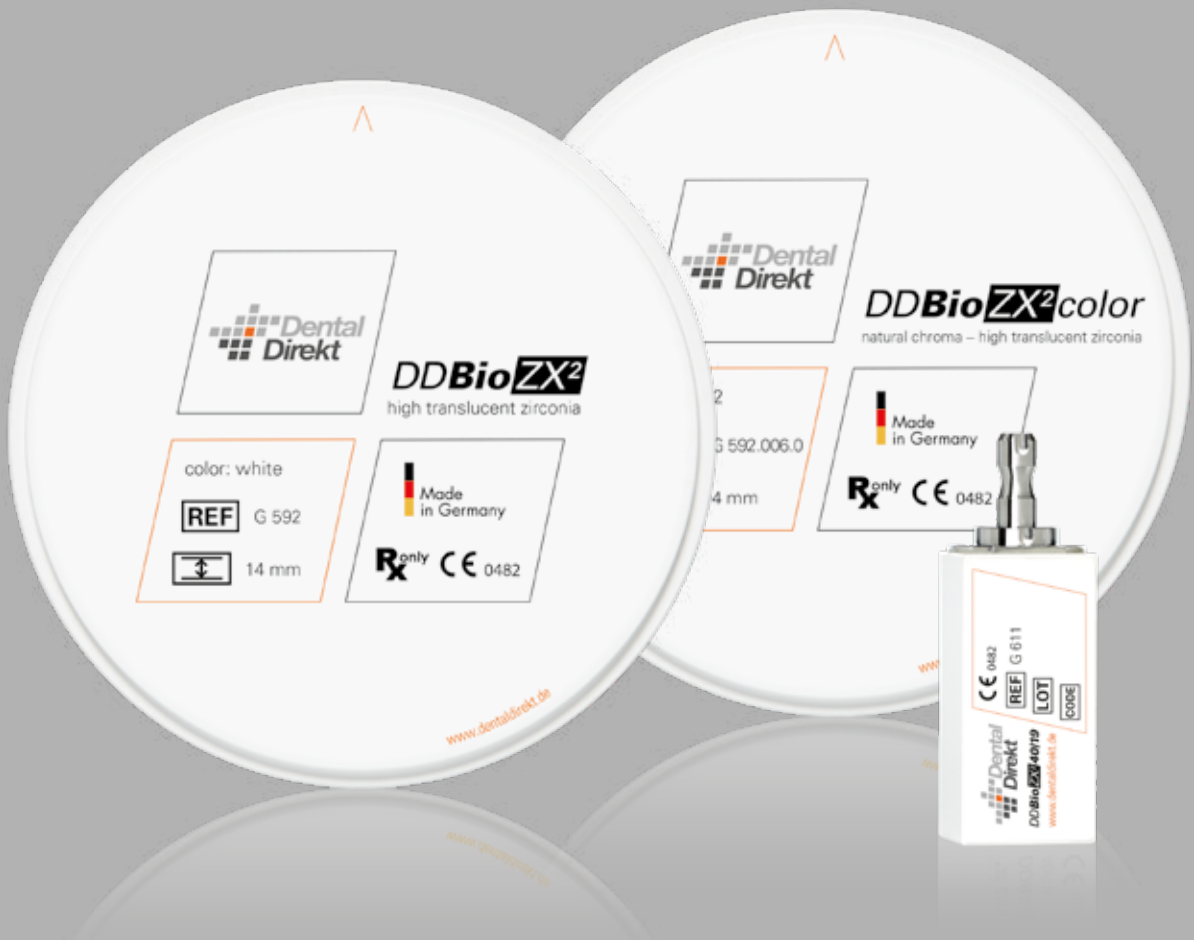

DDBioZX²

high translucent zirconia

High Translucent (3Y-TZP-LA)

WHITE

Benefits

- Further development of the 'classic' zirconia
- High strength
- High fracture toughness

Perfectly suitable for

- Monolithic crowns and bridges (of any span range)
- Especially for monolithic bridges
- Cut back (one layer)
- High aesthetic veneering (of any span range)




Technical data

Type	3Y-TZP-LA
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges with max. 2 connected pontics
Color	White
Aesthetics	High Translucent
Sintering temperature	1.450 °C
Flexural strength	1.150 (±150) MPa*
CAM system	Open systems for Ø 98,5 mm blanks
DD Coloring liquids	DD Pro Shade Z DD Basic Shade DD Art Elements
Aesthetic Finish	DD contrast®
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872

Measure	Order-no.	Price (euro/net)
98,5 x 10 mm	G 590	133,00
98,5 x 12 mm	G 591	153,00
98,5 x 14 mm	G 592	164,00
98,5 x 16 mm	G 593	184,00
98,5 x 18 mm	G 594	205,00
98,5 x 20 mm	G 595	215,00
98,5 x 25 mm	G 596	246,00

 Consult instructions for use.
Not all products are available in every country.

DDBio ZX²19

high translucent zirconia

High Translucent (3Y-TZP-LA)



Benefits

- Further development of the 'classic' zirconia
- High strength
- High fracture toughness

Perfectly suitable for


- Monolithic crowns and bridges (of any span range)
- Especially for monolithic bridges
- Cut back (one layer)
- High aesthetic veneering (of any span range)



Technical data

Type	3Y-TZP-LA
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges with max. 2 connected pontics
Color	White
Aesthetics	High Translucent
Sintering temperature	1.450 °C
Flexural strength	1.150 (±150) MPa*
CAM system	Sirona inLab [®] and inLab [®] MC XL system
DD Coloring liquids	DD Pro Shade Z DD Basic Shade DD Art Elements
Aesthetic Finish	DD contrast[®]
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872

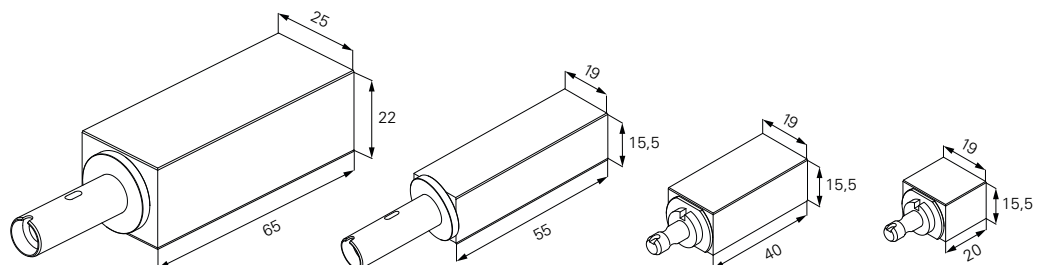

Article	DD Bio ZX ² 65/25	DD Bio ZX ² 55/19 MC XL	DD Bio ZX ² 40/19	DD Bio ZX ² 20/19
Order-no.	G 613	G 612.1	G 611	G 610
Measure in mm	65 x 25 x 22	55 x 19 x 15,5	40 x 19 x 15,5	20 x 19 x 15,5
Item	1	2	10	24
Price (euro/net)	98,00	108,00	273,00	345,00

Consult instructions for use.
Not all products are available in every country.

DDBioZX²color

natural chroma – high translucent zirconia

High Translucent (3Y-TZP-LA)

MONO
CHROME

Benefits

- Further development of the 'classic' zirconia in 5 universal shades
- High strength
- High fracture toughness

Perfectly suitable for

- Monolithic crowns and bridges (of any span range)
- Especially for monolithic bridges
- Cut back (one layer)
- High aesthetic veneering (of any span range)




Technical data

Type	3Y-TZP-LA
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges with max. 2 connected pontics
Color	5 universal colors
Aesthetics	High Translucent
Sintering temperature	1.450 °C
Flexural strength	1.150 (±150) MPa*
CAM system	Open systems for Ø 98,5 mm blanks
Multi Additive Technology®	Reliable reproduction of color via industrial coloring
DD Coloring liquids	DD Art Elements DD Incisal X
Aesthetic Finish	DD contrast®
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872



Measure	light A1 · A2 · B1 B2 · C1	medium A3 · B3 · B4 · C2 C3 · D2 · D3 · D4	intense A3,5 · A4 · C4	low chromatic A1 · A2 · A3 · B2 B3 · D2 · D3	high chromatic A3,5 · A4 · B4	Price (euro/net)
98,5 x 10 mm	G 590.001.0	G 590.002.0	G 590.003.0	G 590.004.0	G 590.005.0	153,00
98,5 x 14 mm	G 592.001.0	G 592.002.0	G 592.003.0	G 592.004.0	G 592.005.0	174,00
98,5 x 18 mm	G 594.001.0	G 594.002.0	G 594.003.0	G 594.004.0	G 594.005.0	215,00
98,5 x 20 mm	G 595.001.0	G 595.002.0	G 595.003.0	G 595.004.0	G 595.005.0	236,00
98,5 x 25 mm	G 596.001.0	G 596.002.0	G 596.003.0	G 596.004.0	G 596.005.0	277,00

 Consult instructions for use.
Not all products are available in every country.

DDBioZX²color

natural chroma – high translucent zirconia according to VITA® shades

High Translucent (3Y-TZP-LA)



Benefits

- Further development of the 'classic' zirconia in 16 VITA® shades
- High strength
- High fracture toughness

Perfectly suitable for

- Monolithic crowns and bridges (of any span range)
- Especially for monolithic bridges
- Cut back (one layer)
- High aesthetic veneering (of any span range)



Technical data

Type	3Y-TZP-LA
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges with max. 2 connected pontics
Color	Precolored acc. to VITA®
Aesthetics	High Translucent
Sintering temperature	1.450 °C
Flexural strength	1.150 (±150) MPa*
CAM system	Open systems for Ø 98,5 mm blanks
Multi Additive Technology®	Reliable reproduction of color via industrial coloring
DD Coloring liquids	DD Art Elements DD Incisal X
Aesthetic Finish	DD contrast®
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872



Measure	A1	A2	A3	A3,5	A4	Price (euro/net)
98,5 x 10 mm	G 590.012.0	G 590.006.0	G 590.007.0	G 590.008.0	G 590.013.0	153,00
98,5 x 14 mm	G 592.012.0	G 592.006.0	G 592.007.0	G 592.008.0	G 592.013.0	174,00
98,5 x 18 mm	G 594.012.0	G 594.006.0	G 594.007.0	G 594.008.0	G 594.013.0	215,00
98,5 x 20 mm	G 595.012.0	G 595.006.0	G 595.007.0	G 595.008.0	G 595.013.0	236,00
98,5 x 25 mm	G 596.012.0	G 596.006.0	G 596.007.0	G 596.008.0	G 596.013.0	277,00

Consult instructions for use.
Not all products are available in every country.



Measure	B1	B2	B3	B4	Price (euro/net)
98,5 x 10 mm	G 590.014.0	G 590.009.0	G 590.015.0	G 590.016.0	153,00
98,5 x 14 mm	G 592.014.0	G 592.009.0	G 592.015.0	G 592.016.0	174,00
98,5 x 18 mm	G 594.014.0	G 594.009.0	G 594.015.0	G 594.016.0	215,00
98,5 x 20 mm	G 595.014.0	G 595.009.0	G 595.015.0	G 595.016.0	236,00
98,5 x 25 mm	G 596.014.0	G 596.009.0	G 596.015.0	G 596.016.0	277,00




Measure	C1	C2	C3	C4	Price (euro/net)
98,5 x 10 mm	G 590.017.0	G 590.010.0	G 590.018.0	G 590.019.0	153,00
98,5 x 14 mm	G 592.017.0	G 592.010.0	G 592.018.0	G 592.019.0	174,00
98,5 x 18 mm	G 594.017.0	G 594.010.0	G 594.018.0	G 594.019.0	215,00
98,5 x 20 mm	G 595.017.0	G 595.010.0	G 595.018.0	G 595.019.0	236,00
98,5 x 25 mm	G 596.017.0	G 596.010.0	G 596.018.0	G 596.019.0	277,00



Measure	D2	D3	D4	Price (euro/net)
98,5 x 10 mm	G 590.020.0	G 590.011.0	G 590.021.0	153,00
98,5 x 14 mm	G 592.020.0	G 592.011.0	G 592.021.0	174,00
98,5 x 18 mm	G 594.020.0	G 594.011.0	G 594.021.0	215,00
98,5 x 20 mm	G 595.020.0	G 595.011.0	G 595.021.0	236,00
98,5 x 25 mm	G 596.020.0	G 596.011.0	G 596.021.0	277,00

DDBioZ

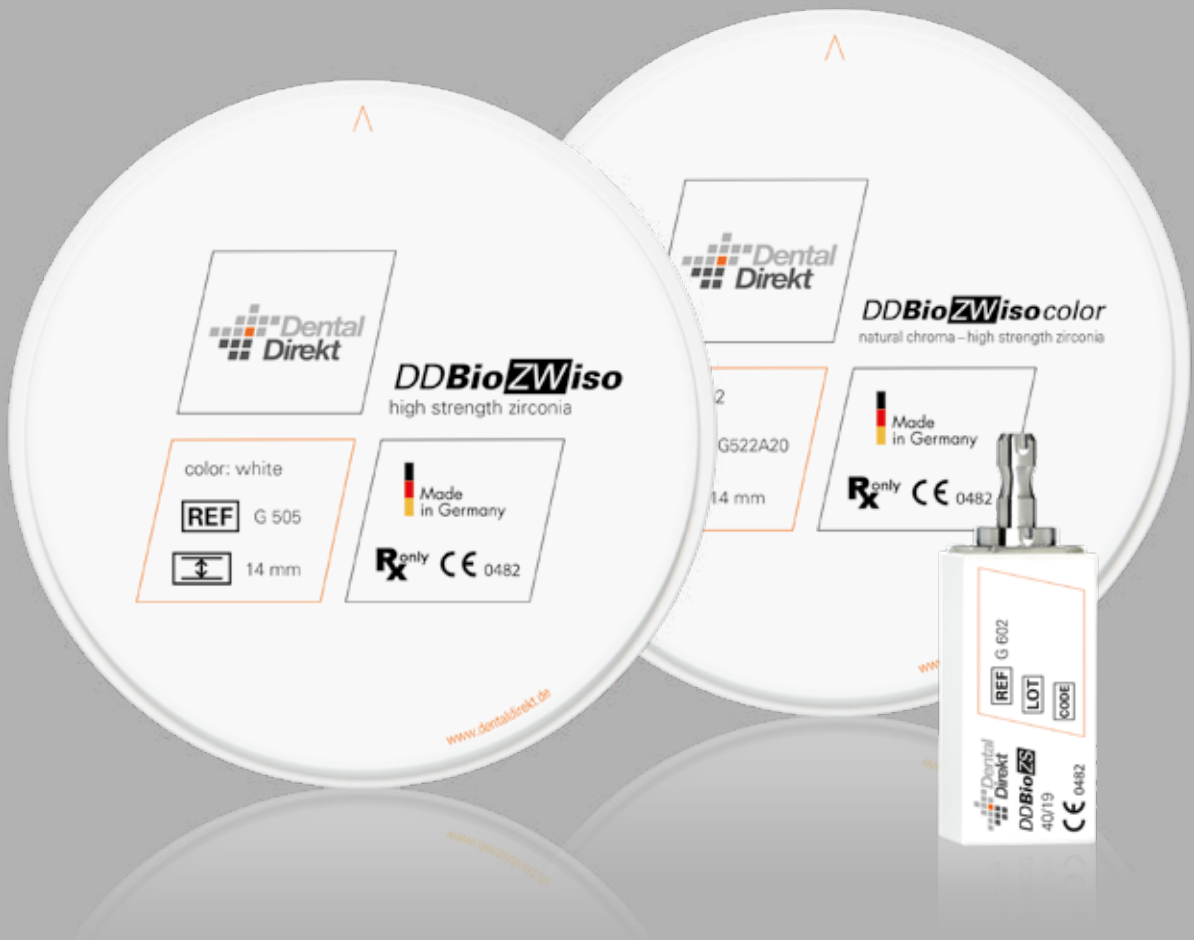

DDBioZWiso

high strength zirconia

High Strength (3Y-TZP-A)

WHITE

Benefits

- The 'classic' among the zirconium oxides
- High strength
- Opaque – ideal for masking
- Very good dyeability

Perfectly suitable for

- Aesthetic veneering of crowns, bridges (of any span range)
- Hybrid abutments




Technical data

Type	3Y-TZP-A
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges with max. 2 connected pontics
Color	White
Aesthetics	Opaque
Sintering temperature	1.450 °C
Flexural strength	1.150 (±200) MPa*
CAM system	Open systems for Ø 98,5 mm blanks
DD Coloring liquids	DD Pro Shade Z DD Basic Shade DD Art Elements
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872

Measure	Order-no.	Price (euro/net)
98,5 x 10 mm	G 513	123,00
98,5 x 12 mm	G 526	143,00
98,5 x 14 mm	G 505	153,00
98,5 x 16 mm	G 527	174,00
98,5 x 18 mm	G 506	195,00
98,5 x 20 mm	G 507	205,00
98,5 x 25 mm	G 553	236,00

 Consult instructions for use.
Not all products are available in every country.

DDBioZS

high strength zirconia

High Strength (3Y-TZP-A)



Benefits

- The 'classic' among the zirconium oxides
- High strength
- Opaque – ideal for masking
- Very good dyeability





Perfectly suitable for

- Aesthetic veneering of crowns, bridges
- Hybrid abutments





Technical data

Type	3Y-TZP-A
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges with max. 2 connected pontics
Color	White
Aesthetics	Opaque
Sintering temperature	1.450 °C
Flexural strength	1.150 (±200) MPa*
CAM system	Sirona inLab® and inLab® MC XL system
DD Coloring liquids	DD Pro Shade Z DD Basic Shade DD Art Elements
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872

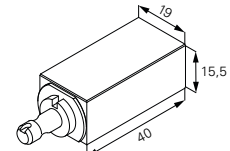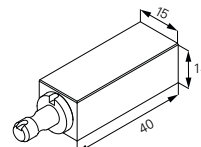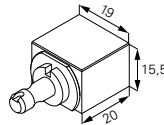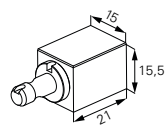

Article	DD Bio ZS 21/15	DD Bio ZS 20/19	DD Bio ZS 40/15	DD Bio ZS 40/19
Order-no.	G 603	G 600	G 601	G 602
Measure in mm	21 x 15 x 15,5	20 x 19 x 15,5	40 x 15 x 14	40 x 19 x 15,5
Item	25	24	10	10
Price (euro/net)	308,00	345,00	258,00	273,00

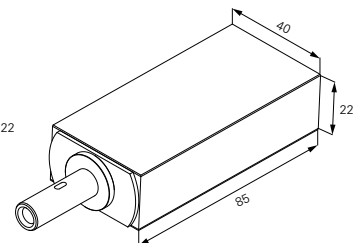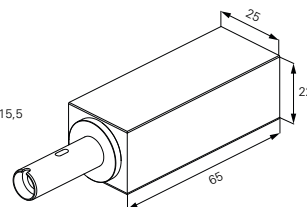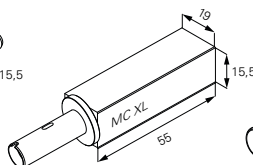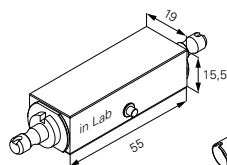

Article	DD Bio ZS 55/19 inLab®	DD Bio ZS 55/19 MC XL	DD Bio ZS 65/25	DD Bio ZS 85/40
Order-no.	G 604	G 604.1	G 605	G 607
Measure in mm	55 x 19 x 15,5	55 x 19 x 15,5	65 x 25 x 22	85 x 40 x 22
Item	2	2	1	1
Price (euro/net)	108,00	108,00	98,00	211,00

Consult instructions for use.
Not all products are available in every country.

DDBioZWiso color

high strength zirconia

High Strength (3Y-TZP-A)

MONO
CHROME

Benefits

- The 'classic' among the zirconium oxides in 5 universal shades
- High strength
- Opaque – ideal for masking

Perfectly suitable for

- Aesthetic veneering of crowns, bridges (of any span range)
- Hybrid abutments



Technical data

Type	3Y-TZP-A
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges s with max. 2 connected pontics
Color	5 universal colors
Aesthetics	Opaque
Sintering temperatur	1.450 °C
Flexural strength	1.150 (±200) MPa*
CAM system	Open systems for Ø 98,5 mm blanks
DD Coloring liquids	DD Art Elements
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872



Measure	500 A1 · A2	800 A3 · B3	1000 C2 · C3	1333 A3,5 · B4	2000 A4	Price (euro/net)
98,5 x 10 mm	G 520.500.0	G 520.800.0	G 520.1000.0	G 520.1333.0	G 520.2000.0	143,00
98,5 x 14 mm	G 522.500.0	G 522.800.0	G 522.1000.0	G 522.1333.0	G 522.2000.0	164,00
98,5 x 18 mm	G 524.500.0	G 524.800.0	G 524.1000.0	G 524.1333.0	G 524.2000.0	205,00
98,5 x 20 mm	G 525.500.0	G 525.800.0	G 525.1000.0	G 525.1333.0	G 525.2000.0	226,00
98,5 x 25 mm	G 526.500.0	G 526.800.0	G 526.1000.0	G 526.1333.0	G 526.2000.0	267,00



Consult instructions for use.

Not all products are available in every country.

DDBioZWisocolor

high strength zirconia
according to VITA® shades

High Strength (3Y-TZP-A)



Benefits

- The 'classic' among the zirconium oxides in 8 VITA®- shades
- High strength
- Opaque – ideal for masking

Perfectly suitable for

- Aesthetic veneering of crowns, bridges (of any span range)
- Hybrid abutments

Technical data	
Type	3Y-TZP-A
Classification	Class IIa according to MDD/MDR
Indication ISO 6872	Type II, class 5 – all dental constructions, bridges with max. 2 connected pontics
Color	Precolored acc. to VITA®
Aesthetics	Opaque
Sintering temperature	1.450 °C
Flexural strength	1.150 (±200) MPa*
CAM system	Open systems for Ø 98,5 mm blanks
DD Coloring liquids	DD Art Elements
Luting system	DD Solid Link

* measured according to DIN EN ISO 6872



Measure	A1	A2	A3	A3,5	Price (euro/net)
98,5 x 10 mm	G520A10	G520A20	G520A30	G520A350	143,00
98,5 x 14 mm	G522A10	G522A20	G522A30	G522A350	164,00
98,5 x 18 mm	G524A10	G524A20	G524A30	G524A350	205,00
98,5 x 20 mm	G525A10	G525A20	G525A30	G525A350	226,00
98,5 x 25 mm	G526A10	G526A20	G526A30	G526A350	267,00



Measure	B1	B2	C1	D2	Price (euro/net)
98,5 x 10 mm	G520B10	G520B20	G520C10	G520D20	143,00
98,5 x 14 mm	G522B10	G522B20	G522C10	G522D20	164,00
98,5 x 18 mm	G524B10	G524B20	G524C10	G524D20	205,00
98,5 x 20 mm	G525B10	G525B20	G525C10	G525D20	226,00
98,5 x 25 mm	G526B10	G526B20	G526C10	G526D20	267,00

Consult instructions for use.
Not all products are available in every country.

Coloring Liquids

02

Coloring Liquid recommendation

exemplified by the shade reference VITA® A2



DDBasicShade
Universal Dentine Liquid



DDBioZ



← more grayish



DDBioZX²



← more grayish



DDcubeONE[®]



← more grayish



DDcubeX²[®]



← more grayish

Color concept:
(shade/color saturation)

intensive chroma, higher gray value

Brightness:


LESS brightness

DD recommendation:

minimum layering (cut back),
full veneering, combination of both techniques

Recommendation



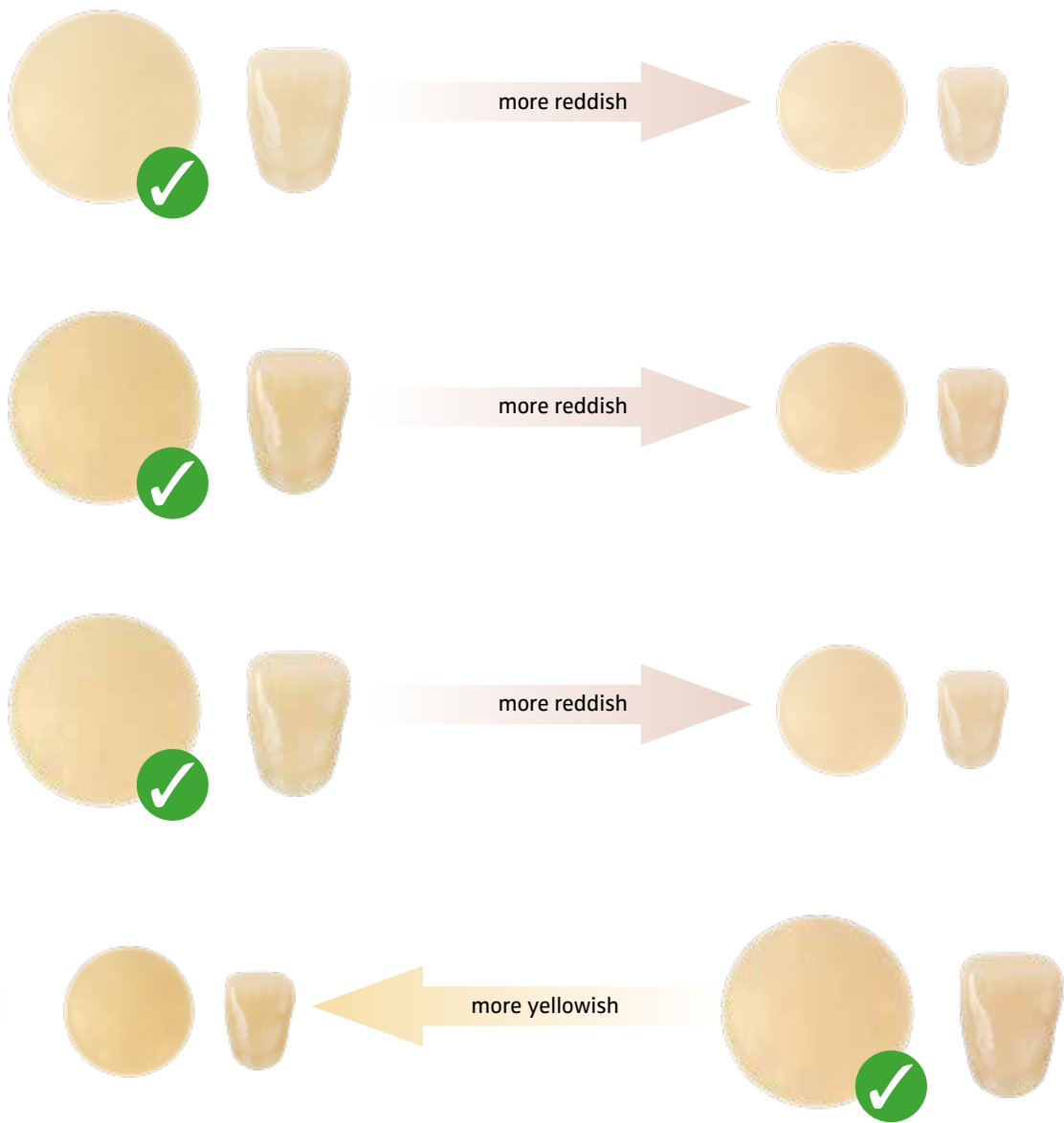

Color concept: (shade/color saturation)

less intensive chroma, warmer shade

Brightness:

MORE brightness

DD recommendation:

Basic coloring with the recommended DD Pro Shade.

Further characterization: For monolithic crowns/bridges made of DD cubeX2® or DD cube ONE® an individual coloring (DD contrast® color and form concept) or minimal layering (cut back). For constructions made of DD Bio ZX2 or DD Bio Z a ceramic veneer.

DDBasicShade

Universal Dentine Liquid

DD Basic Shade is a universal dentine coloring liquid for all non-pre-shaded Dental Direkt zirconia types. DD Basic Shade is used before the sintering process.

Shade guide	VITA® classical (dentine)
Color concept	Intensive chroma – higher gray value
Brightness	Less brightness



16 dentine shades and 3 bleach shades according to VITA®


Color according to VITA®	Order no. 30 ml	Order no. 100 ml
0M1	LBS300M1	LBS1000M1
0M2	LBS300M2	LBS1000M2
0M3	LBS300M3	LBS1000M3
A1	LBS30A1	LBS100A1
A2	LBS30A2	LBS100A2
A3	LBS30A3	LBS100A3
A3,5	LBS30A35	LBS100A35
A4	LBS30A4	LBS100A4
B1	LBS30B1	LBS100B1
B2	LBS30B2	LBS100B2
B3	LBS30B3	LBS100B3
B4	LBS30B4	LBS100B4
C1	LBS30C1	LBS100C1
C2	LBS30C2	LBS100C2
C3	LBS30C3	LBS100C3
C4	LBS30C4	LBS100C4
D2	LBS30D2	LBS100D2
D3	LBS30D3	LBS100D3
D4	LBS30D4	LBS100D4
Price (euro/net)	26,00	66,00

Recommended for

- DD Bio Z (3Y-TZP-A) High Strength
- DD Bio ZX² (3Y-TZP-LA) High Translucent
- DD cube ONE® (4Y-TZP) High Translucent Plus
- DD cubeX²® (5Y-TZP) Super High Translucent

Recommended use

- Color basis for minimal layering (cut back)
- Color basis for full veneering
- Combination of both techniques

 Consult instructions for use.
Not all products are available in every country.

DDProShade Z

Special Dentine Liquid

DD Pro Shade Z is a specially developed dentine coloring liquid for the non-pre-shaded high-strength and high translucent zirconia types. DD Pro Shade Z is used before the sintering process.

Shade guide	VITA® classical (dentine)
Color concept	Less intensive chroma – warmer shade
Brightness	More brightness



16 dentine shades and 3 bleach shades according to VITA®


Color according to VITA®	Order no. 30 ml	Order no. 100 ml
0M1	LPSZ300M1	LPSZ1000M1
0M2	LPSZ300M2	LPSZ1000M2
0M3	LPSZ300M3	LPSZ1000M3
A1	LPSZ30A1	LPSZ100A1
A2	LPSZ30A2	LPSZ100A2
A3	LPSZ30A3	LPSZ100A3
A3,5	LPSZ30A35	LPSZ100A35
A4	LPSZ30A4	LPSZ100A4
B1	LPSZ30B1	LPSZ100B1
B2	LPSZ30B2	LPSZ100B2
B3	LPSZ30B3	LPSZ100B3
B4	LPSZ30B4	LPSZ100B4
C1	LPSZ30C1	LPSZ100C1
C2	LPSZ30C2	LPSZ100C2
C3	LPSZ30C3	LPSZ100C3
C4	LPSZ30C4	LPSZ100C4
D2	LPSZ30D2	LPSZ100D2
D3	LPSZ30D3	LPSZ100D3
D4	LPSZ30D4	LPSZ100D4
Price (euro/net)	26,00	66,00

Recommended for

- DD Bio Z (3Y-TZP-A) High Strength
- DD Bio ZX² (3Y-TZP-LA) High Translucent
- DD cube ONE® (4Y-TZP) High Translucent Plus

Recommended use

- Perfect color base for monolithic crowns and bridges made of DD cube ONE® (DD contrast® color and form concept)
- Color basis for minimal layering (cut back)

 Consult instructions for use.
Not all products are available in every country.

DDProShade^C

Special Dentine Liquid

DD Pro Shade C is a specially developed dentine coloring liquid for the cubic non-pre-shaded zirconia DD cubeX²®. DD Pro Shade C is used before the sintering process.

Shade guide	VITA® classical (dentine)
Color concept	Less intensive chroma – warmer shade
Brightness	More brightness



16 dentine shades and 3 bleach shades according to VITA®

Color according to VITA®	Order no. 30 ml	Order no. 100 ml
0M1	LPSC300M1	LPSC1000M1
0M2	LPSC300M2	LPSC1000M2
0M3	LPSC300M3	LPSC1000M3
A1	LPSC30A1	LPSC100A1
A2	LPSC30A2	LPSC100A2
A3	LPSC30A3	LPSC100A3
A3,5	LPSC30A35	LPSC100A35
A4	LPSC30A4	LPSC100A4
B1	LPSC30B1	LPSC100B1
B2	LPSC30B2	LPSC100B2
B3	LPSC30B3	LPSC100B3
B4	LPSC30B4	LPSC100B4
C1	LPSC30C1	LPSC100C1
C2	LPSC30C2	LPSC100C2
C3	LPSC30C3	LPSC100C3
C4	LPSC30C4	LPSC100C4
D2	LPSC30D2	LPSC100D2
D3	LPSC30D3	LPSC100D3
D4	LPSC30D4	LPSC100D4
Price (euro/net)	26,00	66,00

Recommended for

- DD cubeX²® (5Y-TZP) Super High Translucent

Recommended use

- Perfect color base for monolithic crowns and bridges made of DD cubeX²® (DD contrast® color and form concept)
- Color basis for minimal layering (cut back)



Consult instructions for use.

Not all products are available in every country.

DDArtElements











Recommended for


- DD Bio Z/color (3Y-TZP-A) High Strength
- DD Bio ZX²/color (3Y-TZP-LA) High Translucent
- DD cube ONE®/ML (4Y-TZP) High Translucent Plus
- DD cubeX²®/ML (5Y-TZP) Super High Translucent

To be used with all DD Dentine Liquids.

Effect colors

For individualisation, to shade the dentine colors, or simply for the characterization of e.g. fissures, incisal areas or gingival areas.

		
light pink	purple	blue
LAE30LP	LAE30P	LAE30B
		
yellow	orange	light brown
LAE30Y	LAE30O	LAE30LB
		
dark brown	light grey	dark grey
LAE30DB	LAE30LG	LAE30DG
	30 ml	Price (euro/net)
graphite		
LAE30G		66,00

 Consult instructions for use.
Not all products are available in every country.


Incisal

For occlusal and incisal individualisation.

Color		Order no. / 15 ml
SA1	for painting technique (0M3 · Aa1 · A2 · B1 · B2)	LAE15SA1
SA2	for painting technique (A3 · A3,5 · A4 · B3 · B4)	LAE15SA2
SC1	for painting technique (C1 · C2 · D2)	LAE15SC1
SC2	for painting technique (C3 · C4 · D3 · D4)	LAE15SC2
SO	for painting technique (0M1 · 0M2) for dipping technique (all dentine colors)	LAE15SO
SO Bright	for painting technique (0M1 · 0M2) for dipping technique (all dentine colors) developed specially for DD cube ONE®.	LAE15SOB
Price (euro/net)		17,00

Incisal SO / SO Bright


Optimized for efficient dipping technique.

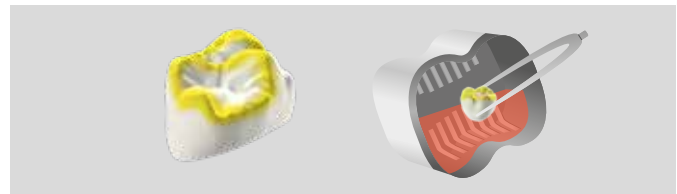

SO Bright produces more intense brightness in the incisal area. Recommended for use with DD cube ONE®.

Incisal SA / SC

Optimized for individual painting technique.

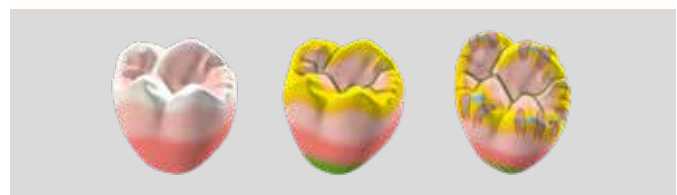

Vario

To brighten the dentine colors.

Article	Order no. / 30 ml
Variation liquid	LAE30V
Price (euro/net)	26,00

DDShadeConcept® sets



Coloring sets	Order no.	Content	Price (euro/net)
DD Basic Shade set	LBS001	16 dentine colors A1 – D4 á 30 ml 5 incisal colors SA1 – SO á 15 ml + 40 applicators	387,00
DD Pro Shade C set (Ideal for DD cubeX2®)	LPSC001	16 dentine colors A1 – D4 á 30 ml 5 incisal colors SA1 – SO á 15 ml + 40 applicators	387,00
DD Pro Shade Z set (Ideal for DD Bio Z, DD Bio ZX2 and DD cube ONE®)	LPSZ001	16 dentine colors A1 – D4 á 30 ml 5 incisal colors SA1 – SO á 15 ml + 40 applicators	387,00
DD Art Elements set	LAE001	5 incisal colors á 15 ml 10 effect colors á 30 ml 1 vario liquid á 30 ml + 40 applicators	608,00



Consult instructions for use.

Not all products are available in every country.

DD Incisal X

Translucency-increasing liquid

The liquid has been specially developed to "pimp" the incisal area of high-strength zirconium oxide restorations and thus ensure even more aesthetic results. DD Incisal X is available in two variants and offers you the greatest flexibility in individualizing your crown and bridge restorations. The application with the enclosed applicators simplifies the handling. By recommended application, the translucent gradient is enhanced.

In short: similarly beautiful to a veneer, but without the great effort.

Benefits

- Two variants for creative working
- Better aesthetics in the incisal area for high-strength, opaque materials
- Aesthetic enhancement without the expense of a veneer

DD Incisal X – X0

Pure translucency-increasing liquid

Add the finishing touch to your restoration. Create contrast between the dentin and the incisal edge or perfect your individualization in combination with DD Art Elements.



Recommended for

- DD cube ONE® ML
- DD Bio ZX² color

Content:


- 1 watch glass
- 10 applicators
- 5 ml DD Incisal X-Liquid (variant X0 or X2) in pipette bottle
- quick guide

DD Incisal X – X2

Translucency-increasing liquid with colouring oxides

Create with minimal effort a real and optical translucency increase through bluish/grayish/violet color effect.

Article	Order No. / 5 ml	Price (euro/net)
DD Incisal X – X0	LIX05X0	66,00
DD Incisal X – X2	LIX05X2	66,00

 Consult instructions for use.
Not all products are available in every country.

Aesthetic Finish

03

DDcontrast®



















color and form concept by *A. Wille*


DD contrast®, the color and texture system from Dental Direkt, gives full monolithic restorations optimum color dynamics, vivid depth effect and various intensities of fluorescence. It creates a wonderful symbiosis with the color and light systems of our monolithic zirconium oxides.



DD contrast® color












- Body, mamelon and effect colors to create characteristic contrasts of light and dark, warm and cold.
- Creation of natural looking tooth effects such as translucency and mamelons.
- With »fissure« the system contains a marking color.

Body colors	Order-no. / 4 g	FG	Price (euro/net)
 dentin A	DDCCDA	F+	47,00
 dentin B	DDCCDB	F+	47,00
 dentin C	DDCCDC	F+	47,00
 dentin D	DDCCDD	F+	47,00
Mamelon colors			
 mamelon maize	DDCCMM	F+	47,00
 mamelon corn	DDCCMC	F+	47,00
 mamelon pink	DDCCMP	F+	47,00
Effect colors			
 white	DDCCWH	F++	47,00
 vanilla	DDCCVA	F+	47,00
 halo	DDCCHA	F+	47,00
 cervical	DDCCCE	F	47,00
 orange	DDCCOR	–	47,00
 traces	DDCCTR	–	47,00
 fissure (marking)	DDCCFI	–	47,00
 light gray	DDCCLG	–	47,00
 dark gray	DDCCDG	–	47,00
 violet	DDCCVI	–	47,00
 blue	DDCCBL	–	47,00

 Consult instructions for use.
Not all products are available in every country.

DD contrast® texture

- Individual modeling of surface structures.
- Vivid depth effect and color intensity on monolithic crowns starting with a minimum layer thickness of 0.1 mm.
- Natural color effect and light dynamics are achieved by a fine feldspar structure.
- Special dentines for increasing the color intensity, of the brightness value and opacity in the labial and buccal area.
- Can be used without additional glaze firing.

Enamel masses	Order-no. / 4 g	FG	Price (euro/net)
 enamel 59	DDCTE59	F+	51,00
 enamel bleach	DDCTEB	F++	51,00
 transpa warm fluor	DDCTWF	F++	51,00
Intensive enamel masses			
 retzius	DDCTRE	F	51,00
 opal blue	DDCTOB	–	51,00
 illusion	DDCTIL	–	51,00
Dentine			
 bleach dentin	DDCTBD	F++	51,00
 chroma dentin A1	DDCTCDA1	F++	51,00
 chroma dentin B4	DDCTCDB4	F+	51,00
Special textures			
 mask	DDCTMA	F	51,00
 caninus	DDCTCA	F	51,00

Fluorescence Grade (FG):   

DD contrast® color, DD contrast® texture and DD contrast® glaze have various intensity of fluorescence. For a better orientation we have designated the level of intensity as follows: F = low fluorescence, F+ = high fluorescence and F++ = very high fluorescence.

DD contrast® glaze

DD contrast® glaze clear (non-fluorescent):

- To seal the surface.
- Glaze with a slight opacity/gray shimmer. Especially advantageous for elderly patients whose teeth are more translucent.

DD contrast® glaze clear fluor (fluorescent):

- To seal the surface.
- Especially in the body area you can achieve the necessary fluorescence – just like with a natural tooth.

Glazes	Order-no. / 4 g	FG	Price (euro/net)
glaze clear	DDCGC	–	43,00
glaze clear fluor	DDCGCF	F++	43,00

DD contrast® liquid

The **DD contrast® standard liquid** can only be used for diluting DD contrast® colors and DD contrast® glaze pastes. With the liquid you can adjust the thickness of the transparent layer.

For the DD contrast® texture pastes, only the **DD contrast® texture liquid** should be used to maintain the special consistency and modelling ability of the material. However, **DD contrast® texture liquid** can also be used with DD contrast® colors and DD contrast® glaze pastes to modify the consistency.

Liquids	Order-no. / 30 ml	Price (euro/net)
standard liquid	DDCSLI	12,50
texture liquid	DDCTLI	12,50


DD contrast® firing paste – the ultimate ceramic firing paste for object fixation

- Stabilizes even difficult ceramic works during the firing process.
- Zero clean-up is required – just a blast of compressed air is sufficient.
- No annoying scraping out or sandblasting after the firing process.
- No drying time is required.



Article	Order-no.	Price (euro/net)
firing paste	DDCFP	22,00

DD contrast® brush	Order-no.	Price (euro/net)
Small – size 3/0	DDCPS	13,00
Medium – size 2	DDCPM	14,00

 Consult instructions for use.
Not all products are available in every country.

DDcontrast® -Set



Article	Order-no.	Content	Price (euro/net)
DD contrast® set	DDCSET	18 x DD contrast® color, à 4 g 11 x DD contrast® texture, à 4 g 1 x DD contrast® glaze clear, 4 g 1 x DD contrast® glaze clear fluor, 4 g 1 x DD contrast® standard liquid, 30 ml 1 x DD contrast® texture liquid, 30 ml 1 x DD contrast® brush small, size 3/0 1 x DD contrast® brush medium, size 2	1.418,00



Consult instructions for use.

Not all products are available in every country.

 **panther**
Das Premium-Finish

Enables polishing of finalized zirconia restorations. In addition to the polishers for initial smoothing (purple), a second variant of polishers (white) is available for gloss, as well as goat hair brushes and a diamond polishing paste.



panther starter kit

Order-no.	TCPSET104/01
Item	1
Price (euro/net)	216,00



TCP M	panther Flame 055	panther Lense 260	panther Wheel 170	panther Lense 150
Order-no.	TCP100104MSW4Z	TCP140104M	TCP158104MSW4Z	TCP541104MSW4Z
Measure (mm)	Ø 5,5 x 18	Ø 26 x 2	Ø 17 x 2,4	Ø 15 x 1
Shaft	104-HP	104-HP	104-HP	104-HP
🔄 Optimal / maximum (min. -1)	12.000 / 20.000	8.000 / 20.000	8.000 / 20.000	8.000 / 20.000
Item	3 pieces	1 piece	3 pieces	3 pieces
Price (euro/net)	25,00	36,00	66,00	72,00



TCP XF	panther Flame 055	panther Lense 260	panther Wheel 170	panther Lense 150
Order-no.	TCP100104XFSW4Z	TCP140104XF	TCP158104XFSW4Z	TCP541104XFSW4Z
Measure (mm)	Ø 5,5 x 18	Ø 26 x 2	Ø 17 x 2,4	Ø 15 x 1
Shaft	104-HP	104-HP	104-HP	104-HP
🔄 Optimal / maximum (min. -1)	12.000 / 20.000	8.000 / 20.000	8.000 / 20.000	8.000 / 20.000
Item	3 pieces	1 piece	3 pieces	3 pieces
Price (euro/net)	25,00	36,00	66,00	72,00



Accessories	panther Diamond polishing paste	Goat's hair brush 17 mm	Goat's hair brush 19 mm
Order-no.	TCPPASTE01/20	TCP417104GR16	TCP419104GR16
Measure (mm)	-	Ø 17	Ø 19
Shaft	-	104-HP	104-HP
🔄 Optimal / maximum (min. -1)	-	12.000 / 15.000	12.000 / 15.000
Item	20 g	10 pieces	10 pieces
Price (euro/net)	78,00	13,00	13,00

High Performance Polymers

04

DDBioSplint^{PHI}

high impact thermoplastic


Polymethylmethacrylate (PMMA)

PMMA for occlusal splints

DD Bio Splint P HI are transparent thermoplastic blanks based on polymethylmethacrylate (PMMA) and suitable for the fabrication of bite splints. In the CAD/CAM process, you achieve reproducible results with little post-processing effort.

Benefits


- Biocompatible
- High stability due to high impact modification
- Very good polishing properties
- Possible to bond to conventional polymers
- Bisphenol A free

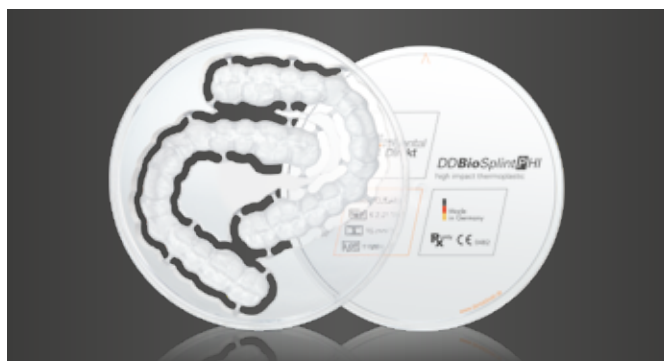

Technical data

Material	100% Polymethylmethacrylate (PMMA)
Classification	Class IIa according to MDD/MDR
Indication	Splints, therapeutic splints, bite regulators (indicated for the long-term use in the oral cavity)
Color	clear
Flexural strength ISO 20795-2	64 MPa
Flexural modulus ISO 20795-2	≥ 1950 MPa
Fracture toughness ISO 20795-2	2,0 MPa m ^{1/2}
Water absorption ISO 20795-2	< 24 µg/mm ³
Solubility ISO 20795-2	< 0,13 µg/mm ³
Residual monomer ISO 20795-2	< 0,7 %
CAM system	open systems for Ø 98,5 mm blanks

Measure	Order-no.	Price (euro/net)
98,5 x 15 mm	K40115	35,00
98,5 x 20 mm	K40120	42,00
98,5 x 25 mm	K40125	47,00
98,5 x 30 mm	K40130	58,00

 Consult instructions for use.
Not all products are available in every country.

DDBioSplintC

innovative polycarbonate

Polycarbonate (PC)

Polycarbonate for occlusal splints

DD Bio Splint C blanks consist of a biocompatible polycarbonate for use in medical devices. In the CAD/CAM process, you achieve reproducible results with little post-processing effort.

Benefits

- Biocompatible
- Low level of flexibility
- Extremely high stability
- Very good polishing properties
- Possible to bond to conventional polymers



Technical data

Material	100% Polycarbonate (PC)
Classification	Class I according to MDD/MDR
Indication	Splints, therapeutic splints, bite regulators (indicated for the short-term use in the oral cavity)
Color	clear/slightly bluish
Flexural strength ISO 20795-1	≥ 98 MPa
Flexural modulus ISO 20795-1	> 2.300 MPa
Water absorption ISO 20795-1	3,6 µg/mm ³
Solubility ISO 20795-1	0,2 µg/mm ³
residual monomer content ISO 20795-1	≤ 0,0015 %
CAM system	open systems for Ø 98,5 mm blanks

Measure	Order-no.	Price (euro/net)
98,5 x 15 mm	K 2.01.15	47,00
98,5 x 20 mm	K 2.01.20	53,00



Consult instructions for use.

Not all products are available in every country.

DDBioSplintFLEX

flexible thermoplastic


Copolyester (PCT-G)

Copolyester for bite splints

The DD Bio Splint FLEX blanks are made of bio-compatible copolyester for use in medical devices. The material is characterised by high flexibility and simultaneously high impact strength. In the CAD/CAM process, you achieve reproducible results with little post-processing effort.

Benefits

- Biocompatible
- High level of flexibility
- Good polishing properties
- Possible to bond to conventional polymers (with bonder)
- Bisphenol A free

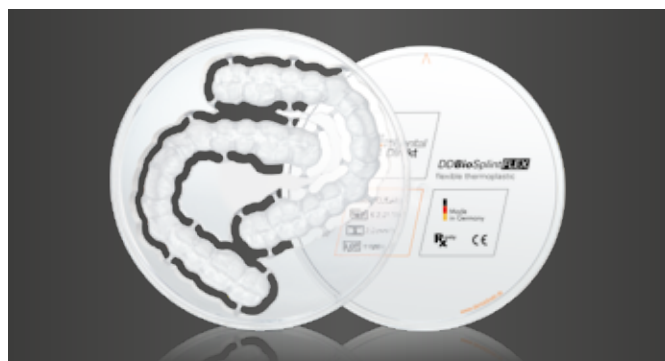

Technical data

Material	100% Copolyester (PCT-G)
Classification	Class I according to MDD
Indication	Flexible splints, therapeutic splints, bite regulators (Indicated for short-term use in the oral cavity)
Color	clear
Flexural strength ISO 20795-1	≥ 58 MPa
Flexural modulus ISO 20795-1	≥ 1.370 MPa
Fracture toughness ISO 20795-1	≥ 3 MPa m ^{1/2}
Water absorption ISO 20795-1	6 µg/mm ³
Solubility ISO 20795-1	< 0,25 µg/mm ³
CAM system	open systems for Ø 98,5 mm blanks

Measure	Order-no.	Price (euro/net)
98,5 x 20 mm	K60120	63,00



Consult instructions for use.

Not all products are available in every country.

DDproviPHI

high impact thermoplastic

Polymethylmethacrylate (PMMA)



PMMA for temporary restorations

DD provi P HI is a monochromatic colored milling blank and based on the raw material of our proven DD Bio Splint P HI. The material is perfectly balanced in terms of translucency and opacity and adapted to the VITA® Classic shades.

Benefits

- Biocompatible
- High stability due to high impact modification
- Very good polishing properties
- Customizable
- Possible to bond to conventional polymers



Technical data

Material	Polymethylmethacrylate (PMMA)
Classification	Class IIa according to MDD/MDR
Indication	Temporary crowns and bridges (indicated for the long-term use in the oral cavity)
Flexural strength ISO 20795-1	65 MPa
Flexural modulus ISO 20795-1	2.030 MPa
Fracture toughness ISO 20795-1	2,6 MPa m ^{1/2}
Water absorption ISO 20795-1	< 24 µg/mm ³
Solubility ISO 20795-1	< 0,3 µg/mm ³
Residual monomer content ISO 20795-1	< 0,5 %
CAM system	open systems for Ø 98,5 mm blanks

Chemical Composition [Weight %]

Polymethylmethacrylate (PMMA)	> 99
Dyes	< 1



Measure	Bleach	A1	A2	A3	A3,5	B2	C2	D3	Price (euro/net)
98,5 x 20 mm	K71720	K70120	K70220	K70320	K70420	K70720	K71120	K71520	42,00

Consult instructions for use.
Not all products are available in every country.

DDbase^{PHI}

high impact thermoplastic

Polymethylmethacrylate (PMMA)



PMMA for digital denture bases

DD base P HI are gingiva-colored milling blanks based on the raw material of our proven DD Bio Splint P HI. The material is ideally suited for the production of digital denture bases and indicated for the long-term use in the oral cavity.


Benefits

- Biocompatible
- High stability due to high impact modification
- Very good polishing properties
- Very good bond to denture teeth
- Available in four colors



Technical data	
Material	Polymethylmethacrylate (PMMA)
Classification	Class IIa according to MDD/MDR
Indication	Denture bases for removable dentures (indicated for the long-term use in the oral cavity)
Flexural strength ISO 20795-1	65 MPa
Flexural modulus ISO 20795-1	2.030 MPa
Fracture toughness ISO 20795-1	2,6 $\mu\text{g}/\text{mm}^{1/2}$
Water absorption ISO 20795-1	< 24 $\mu\text{g}/\text{mm}^3$
Solubility ISO 20795-1	< 0,3 $\mu\text{g}/\text{mm}^3$
residual monomer content ISO 20795-1	< 0,5 %
CAM system	open systems for \varnothing 98,5 mm blanks
Chemical Composition [Weight %]	
Polymethylmethacrylate (PMMA)	> 99
Dyes	< 1

High Performance Polymers

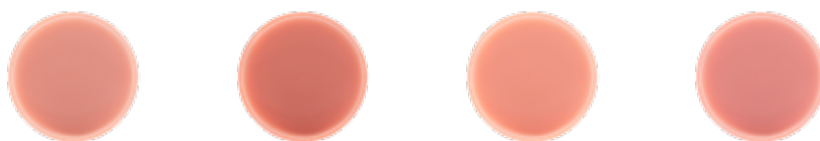

Measure	light pink	pink	sun	rouge	Price (euro/net)
98,5 x 30 mm	K50230	K50130	K50330	K50430	79,00

Consult instructions for use.
Not all products are available in every country.

DDpolyX^{ML}

multilayer – natural colored polymers

Polymethylmethacrylate (PMMA)



Acrylic Polymer

DD poly X ML blanks consist of multilayer acrylic polymer based on polymethylmethacrylate (PMMA).

Benefits

- Biocompatible – medical polymer material
- 5 color layers for a natural aesthetic effect
- Excellent hardness and abrasion resistance
- Good polishing properties
- Customizable with composite



Technical data

Material	Polymethylmethacrylate (PMMA)
Classification	Class IIa according to MDD/MDR
Indication	Implant-supported temporary restorations, temporary crowns and bridges with up to two connected pontics in the front and side tooth area (maximum recommended usage without check-up by dentist is 12 months)
Flexural strength ISO 10477	90 MPa
Flexural modulus ISO 20795-1	2.300 MPa
Water absorption ISO 10477	25 µg/mm ³
Solubility ISO 10477	0,0 µg/mm ³
residual monomer content ISO 20795-1	< 2 %
CAM system	open systems for Ø 98,5 mm blanks
Chemical Composition [Weight %]	
Polymethylmethacrylate (PMMA)	> 99
Coloring pigments	< 1



Measure	B0	A1	A2	A3	A3.5	A4	Price (euro/net)
98,5 x 16 mm	KML0416	KML0116	KML0216	KML0316	KML0716	KML0816	79,00
98,5 x 20 mm	KML0420	KML0120	KML0220	KML0320	KML0720	KML0820	84,00



Consult instructions for use.

Not all products are available in every country.



Measure	B1	B2	B3	B4	Price (euro/net)
98,5 x 16 mm	KML0916	KML1016	KML0516	KML1116	79,00
98,5 x 20 mm	KML0920	KML1020	KML0520	KML1120	84,00



Measure	C1	C2	C3	C4	Price (euro/net)
98,5 x 16 mm	KML1216	KML0616	KML1316	KML1416	79,00
98,5 x 20 mm	KML1220	KML0620	KML1320	KML1420	84,00



Measure	D2	D3	D4	Price (euro/net)
98,5 x 16 mm	KML1516	KML1616	KML1716	79,00
98,5 x 20 mm	KML1520	KML1620	KML1720	84,00

DDpeekMED

the polymer steel

Polyether ether ketone



Medical PEEK

DD peek MED are biocompatible milling blanks based on the high-performance polymer polyetheretherketone.

Benefits

- For allergy patients suitable
- Free from monomer
- Extremely loadable
- Several colours for different applications

Technical data	
Material	Polyether ether ketone (PEEK)
Classification	Class IIa according to MDD/MDR
Indication	For the manufacture of fixed or removable dental framework constructions for long-term use. e.g. abutments and cervical gingival formers, bar constructions, tertiary constructions, telescope crowns, transversal links (indicated for the long-term use in the oral cavity)
Flexural strength* ISO 20795-1	≥ 155 MPa
Flexural modulus* ISO 20795-1	≥ 3,8 GPa
Density* ISO 1183	1,5 g/cm ³
Water absorption* ISO 20795-1	0,0 µg/mm ³
Solubility* ISO 20795-1	0,0 µg/mm ³
CAM system	open systems for Ø 98,5 mm blanks

*Value determination based on DD peek MED natural

Chemical Composition [Weight %]				
	natural	white	dentine	gingiva
Polyether ether ketone (PEEK)	≥ 99,99	≥ 80	≥ 80	> 90
Titanium oxide	–	< 20	< 20	< 5
Further Additives	< 0,01	< 0,1	< 1	< 1



Measure	white	dentine	natural	gingiva	Price (euro/net)
98,5 x 12 mm	K 3.01.12	K 3.02.12	K 3.03.12	K 3.04.12	190,00
98,5 x 14 mm	K 3.01.14	K 3.02.14	K 3.03.14	K 3.04.14	210,00
98,5 x 16 mm	K 3.01.16	K 3.02.16	K 3.03.16	K 3.04.16	230,00
98,5 x 18 mm	K 3.01.18	K 3.02.18	K 3.03.18	K 3.04.18	250,00
98,5 x 20 mm	K 3.01.20	K 3.02.20	K 3.03.20	K 3.04.20	270,00
98,5 x 23 mm	K 3.01.23	K 3.02.23	K 3.03.23	K 3.04.23	290,00
98,5 x 30 mm	K 3.01.30	K 3.02.30	K 3.03.30	K 3.04.30	350,00



Consult instructions for use.

Not all products are available in every country.

More Materials

05

DDBioCW


Milling blank based on CoCr

The DD Bio CW blanks provide a high level of working safety and guarantees reliably good quality. Outstanding aesthetic results can be achieved by a veneering with a material-specific veneering ceramic. The adherence oxide formers provide an optimal bond between metal and ceramic.

The DD Bio CW blank has been treated with a special heat treatment, which makes the alloy easy to machine and gives it outstanding properties due to a uniform microstructure.

Benefits

- Corrosion resistant
- Free from nickel, cadmium and beryllium
- Optimum metal-ceramic bond

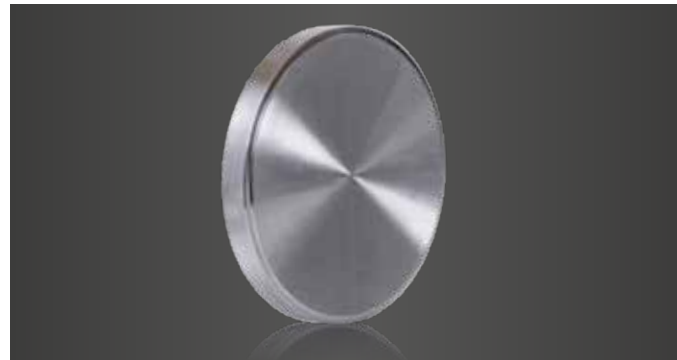

Technical data

Density	8,8 g/cm ³
E modulus	250 GPa
CTE (25 – 500°C)	14,4 10 ⁻⁶ K ⁻¹
Vickers hardness	285 HV10
Elongation at break	9 %
CAM system	open systems for Ø 98,5 mm blanks
Luting system	DD Solid Link

Chemical composition [Weight %]

Co	63
Cr	29
Mo	6
Nb, Si, Mn, Fe	< 1

Measure	Order-no.	Price (euro/net)
98,5 x 8 mm	G 900	175,00
98,5 x 10 mm	G 901	195,00
98,5 x 12 mm	G 902	215,00
98,5 x 13,5 mm	G 903	235,00
98,5 x 15 mm	G 904	245,00
98,5 x 18 mm	G 905	285,00



Consult instructions for use.

Not all products are available in every country.

DDcamWAX

Wax for casting technique

Synthetic wax blanks utilised for the production of combustible wax frameworks that leave no residue. DD CAM wax blanks are distinguished by excellent milling properties and form a fine chip during the machining that does not clog milling cutters (non-medical product).



Specification

CAM system open systems for Ø 98,5 mm blanks

Measure	Order-no.	Color	Price (euro/net)
98,5 x 20 mm	G 529	gray	32,00

DDcalibrationdisc


The DD calibration disc made of polyurethane is used to calibrate your milling machine. After milling a special calibration body, you measure them and enter these values into the program. The milling machine will guide you through the calibration process step by step (non-medical product).



Specification

CAM system open systems for Ø 98,5 mm blanks

Measure	Order-no.	Color	Price (euro/net)
98,5 x 15 mm	G 944	brown	22,00

 Consult instructions for use.
Not all products are available in every country.

Accessories

06

DDShadeGuide

Make your own DD Shade Guide

With the DD Shade Guide, an individual and practice-oriented shade guide is prepared, which is tailored to the zirconium oxides used and the workflow in the dental laboratory. In this way, the tooth shade is determined directly with the respective restoration material. This ensures a much better expectation of the shade effect and improves the reproducibility of the same. This minimizes the risk of a prosthetic failure due to incorrect shade taking.



Example: DD Shade Guide (set of 8 colors) with DD Bio ZW iso color

Benefits

- Shade guide made of the zirconia material of your choice
- 100 % reproducibility
- No transfer errors during color matching due to your own, individual color key

Set content



Sticks and fixation clips
(number according to colors)




DD Shade Guide holder
(Set of 8 or 16)



Stl. data record on USB stick

DD Shade Guide	Order-no.	Price (euro/net)
Set of 8	Z 017	46,00
Set of 16	Z 016	67,00

 Consult instructions for use.
Not all products are available in every country.

DDSolidLink

Universal dual-curing composite based luting for intra and extraoral use

The DD Solid Link is a permanent luting solution for all common dental materials. Due to the reliable and precise luting is the fixation guaranteed. DD Solid Link forms the link between the sub- and suprastructures.

Benefits

- Permanent luting for the common dental materials
- Universal solution, for intra- and extraoral usage
- Perfect masking properties for discolored or titanium bases with highly opaque shades
- Working time approx. 2 min.
- Setting time extraoral (23 °C) approx. 7 min.
- Setting time intraoral (37 °C) approx. 4 min.

Three shades for the perfect color matching

The DD Solid Link is available in two highly opaque (white/high opaque and dentine/high opaque) colors and additionally in a translucent shade.

Content

- 10 g DD Solid Link
- 10 Mixing tips type 9
- Instructions for use



Technical data

Classification	Class IIa according to MDD
Indication	Permanent extra- and intraoral luting of suprastructures (incl. mesostructures) made of metal, ceramics, zirconium dioxide and polymers (composite, hybrid-ceramics or PMMA) <ul style="list-style-type: none"> • Prefabricated adhesive/titanium bases (only extraorally) • Customised titanium and zirconium dioxide abutments (extra- and intraorally) • Tooth structure
Chemical composition (in descending order according to content)	Barium aluminum borosilicate glass, Titanium dioxide, HEDMA, BisGMA, fluorosilicate glass, pyrogenic silicon dioxide, initiators, stabilizers, color pigments
3-Point flexural strength* ISO 4049	119 MPa
Shear Bond strength on zirconium dioxide** ISO/TS 29022	24.4 MPa
Shear Bond strength on titanium** ISO/TS 29022	25.1 MPa
Shear Bond strength on lithium disilicate** ISO/TS 29022	21.7 MPa

* after 24 h storage at 23° C

** with DD Uni Bond

Luting composite

Article	Order-no. / 10 g	Price (euro/net)
DD Solid Link translucent	DDSLCT	108,00
DD Solid Link dentine high opaque	DDSLCD	108,00
DD Solid Link white high opaque	DDSLCW	108,00




Consult instructions for use.

Not all products are available in every country.

DDSolidLinkTry-In

For the precise shade matching

In order to achieve a perfect shade matching there are different Try-In Pastes to use in advance. It gives the possibility for highly aesthetic restorations to check the shades before the permanent fixation. Afterwards the pastes are easy to remove simply just with water.

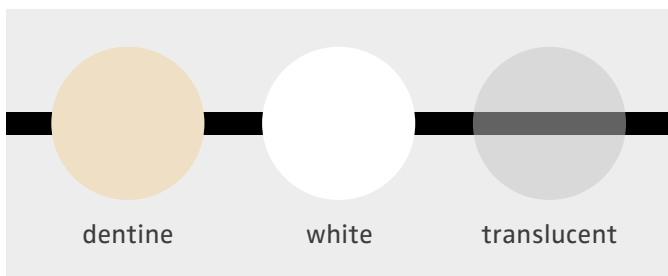

Technical data


Classification	Class I according to MDD
Indication	Try-In Paste for shade matching With Try-In pastes – precise shade simulation for hybrid abutments, abutment crowns and core treatments

Content

- 2,5 g DD Solid Link Try-In Paste
- 10 Mixing tips type 41
- Instructions for use

Try-In

Article	Order-no. / 2,5g	Price (euro/net)
DD Solid Link Try-In translucent	DDSLTIT	38,00
DD Solid Link Try-In dentine high opaque	DDSLTID	38,00
DD Solid Link Try-In white high opaque	DDSLTIW	38,00

 Consult instructions for use.
Not all products are available in every country.

DDUniBond

Bonding agent for ceramics, metal and composite

The one-component DD Uni Bond coupling agent provide a durable adhesive bonding between indirect restorations and DD Solid Link.

Benefits

- Creates optimal conditions for the permanent link
- Easy to use – no mixing necessary

Content


- 5 ml DD Uni Bond
- Instructions for use



Technical data

Classification	Class IIa according to MDD
Indication	DD Uni Bond is used as a coupling agent to achieve an optimal bond strength between indirect restorations made of either silicate ceramics, aluminium oxide, zirconium dioxide, metal or composite, as well as glass fibre-reinforced composite posts, on the one side and methacrylate-based luting composites on the other side.
Chemical composition	Organic acid, 3-methacryloxy-propyltrimethoxysilane and acetone

Article	Order-no. / 5 ml	Price (euro/net)
DD Uni Bond	DDSLUB	77,00

 Consult instructions for use.
Not all products are available in every country.

DDUniBondD

Bonding agent for the tooth structure

The DD Uni Bond D is a dual-curing, self-etching nano-composites reinforced bond for creating a durable bonding between the tooth structure and DD Solid Link.

The DD Uni Bond D consists of 2 components.

Benefits

- Chemical bonding for areas where is light-curing not possible
- Eliminates the need for a separate etching step yet still achieves the same adhesive properties as total-etch bonds
- Tolerates residual moisture
- Durable marginal fitting

Content


- 4 ml Liquid 1
- 4 ml Liquid 2
- 50 single use brush / applicators
- 2 Mixing trays
- Instructions for use



Technical data

Classification	Class IIa according to MDD
Indication	Direct self-curing or dual-curing composite restorations and core build-ups Direct light-curing composite/ compomer/Ormocer®-based restorations Fixing of root pins/posts with dual-curing or self-curing composite cements Indirect restorations; when using dual-curing/self-curing composite cements for bonding of inlays, onlays, crowns and bridges
Chemical composition	Organic acids, Bis-GMA, HEMA, TMPTMA, campherchinon, amines (DABE), BHT, catalysts, fluoride and ethanol.

Article	Order-no. / 2 x 4 ml	Price (euro/net)
DD Uni Bond D	DDSLUBD	123,00

 Consult instructions for use.
Not all products are available in every country.

DDphoeniX

reviving cleaning powder

Despite careful consideration of all important aspects of the laboratory process, discoloration and contamination of the objects may occur. A frequent reason for this are oxides from the heating elements and their deposition in the furnace chamber.

Due to its special granulation, the DD phoeniX cleaning powder has a very reactive surface, which attracts the oxides in the furnace chamber and binds them permanently. Thus supports efficiently the cleaning process.

Benefits


- No molybdenum-caused green discoloration of the sintered objects
- Binds the contaminations from the furnace atmosphere and supports the cleaning
- Regular use increases the color consistency and reproducibility of your results



Important note

Among the high-temperature furnaces for dental use mainly furnaces with heating elements made of molybdenum disilicide (MoSi_2) are affected by oxidic precipitation and discoloration. The sintering furnaces available from Dental Direkt heat with siliciumcarbide heating elements (SiC). Discoloration or contamination are not caused by SiC heating elements.

Article	Order-no.	Price (euro/net)
DD phoeniX bottle à 200 g	Z 010	74,00

 Consult instructions for use.
Not all products are available in every country.

Dental Direkt GmbH
Industriezentrum 106 – 108 | 32139 Spenge | Germany
Tel: +49 5225 86319-0 | E-Mail: info@dentaldirekt.de | www.dentaldirekt.de/en

Rev. 20_2024/02

Art. no.: PM017ENG

CLOSE TO YOU

